



Hyperbaric Oxygen Chamber

(Enhanced 2ATA Soft Body laying chamber for single person)

Model: uDR L3

About Hyperbaric Oxygen Therapy (HBOT):

Hyperbaric oxygen therapy (HBOT) is well known for treating scuba and deep-sea divers affected by the rapid change in pressure around them. But did you know that HBOT is also used to treat a variety of other health problems, including carbon monoxide poisoning and diabetic foot ulcers.

Hyperbaric oxygen therapy treats chronic wounds and other medical conditions. Room air contains 21% oxygen, while hyperbaric oxygen therapy provides you with air that contains 100% oxygen. A pressurized chamber delivers oxygen that's two to three times higher than normal air pressure, which helps your lungs gather and absorb more oxygen. Oxygen helps the tissues in your body heal and withstand infections.

Healthcare providers use hyperbaric oxygen therapy alone or with other medical treatments, such as medications or surgery. Most people receive hyperbaric oxygen therapy on an outpatient basis for one to two hours each session. Depending on the type of treatment, you may sit or lie down during the procedure.

Other names for hyperbaric oxygen therapy include HBO therapy or HBOT.


Features:

- made with enhanced material and special design, it can bear 2ATA pressure. The first 2ATA soft chamber in the world.
- Super Large space inside, it can enjoy together with your friends
- you can lay or sit inside, relax better while inhaling oxygen
- The cabin is made of environmentally friendly degradable TPU, laser cutting, seamless heat sealing, with observation window
- Using YKK zipper, which can be opened in both directions, sealed with double-layer flexible adaptive pressure difference space isolation sealing method




Specification:



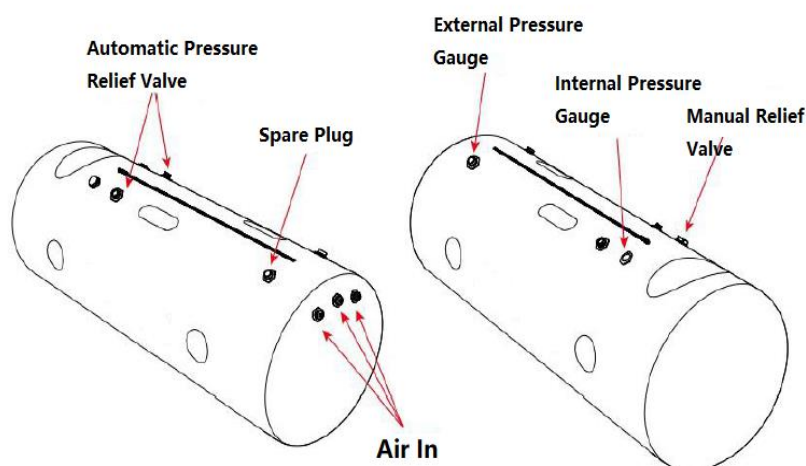


About Chamber:

Model Name	uDR L3
Control System	Touch screen control(10 inch)
Material	TPU
Color	White, blue, or customized
Size	Φ80cm*200cm or customized
Configuration	As the list below
Diffuse oxygen concentration	<45%
Working pressure in the cabin	100-200KPa adjustable
Working noise	<30db
Temperature in the cabin	Ambient temperature +3°C (without air conditioner)
Safety Facilities	Manual safety valve, automatic safety valve
Capacity	One adjustable recliner
Application	Home,Hospital,oxygen bar,oxygen therapy center
Function	hyperbaric oxygen therapy

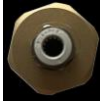

Accessories


Parts Name	Standard Quantity	Picture
Manual relief valve	1PCS	
Automatic pressure relief valve	2PCS	
Air intake	3PCS	

Plug	1PCS	
Silencer	1PCS	
inflatable column Base	2PCS	
Internal pressure gauge	1PCS(optional)	
External pressure gauge	1PCS(optional)	
External bracket	1PCS(optional)	


About Oxygen Supply System:

Model	uMR O8
Size	H768*L420*W400mm
Control System	Touch screen control(10 inch)
Power Supply	AC 100V-240V 50/60Hz
Power	750W
Oxygen Pipe Diameter	8 mm
Air Pipe Diameter	12 mm
Oxygen flow	15L/min
Max airflow	120 L/min
Max Outlet pressure	200KPA(2ATA)
Oxygen Purity	96%±3%(when the chamber is in a fully sealed)
Oxygen system	air filter (PSA)
Compressor	Oil-free compressor air delivery system
Noise	≤45db

UI Interface:

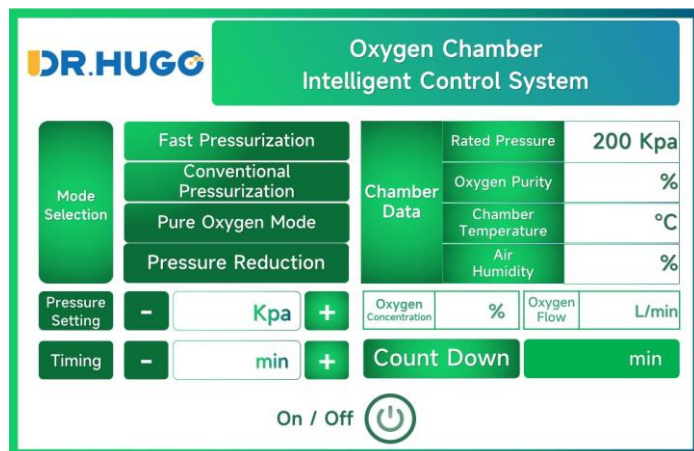

Machine:

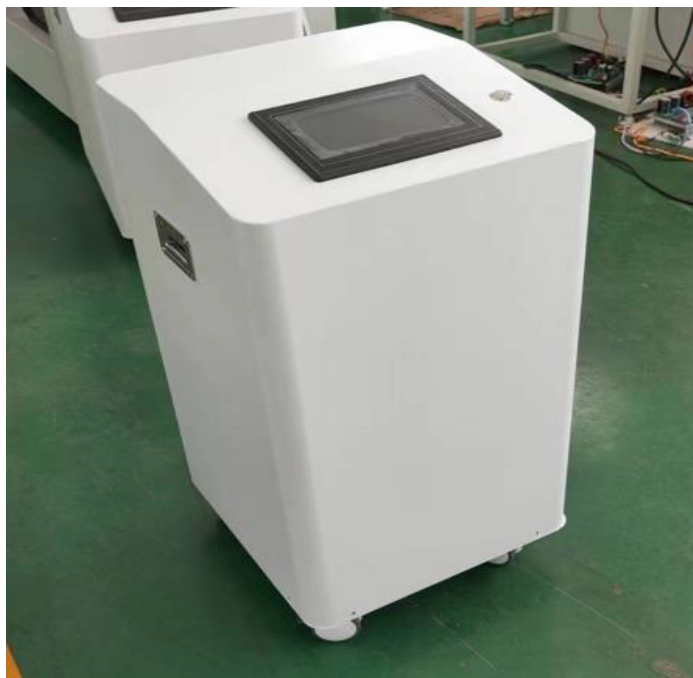

Remot control panel:



Configuration List:

Index	Specification	Unit	QTY
Soft Body Chamber	Φ80cm*200cm or Customize	set	1
Oxygen Supply System	Model: uMR 08 or upgrade to other model(Optional)	set	1
Automatic Safety Relief Valve	3.5Psi/4.5Psi/Customize	set	2
Emergency Relief Valve	standard	set	1
In-cabin controller (touch screen version)	Standard: integrated pressure, oxygen concentration, temperature And can be customize	set	1
Communication System	Optional: Internal and external, Dual wireless, Chargeable	set	1
In-cabin Power Supply System	Optional: 12/24V(double voltage) 600W	set	1
Air Conditioner	Optional : water cool	set	1
Air Cooler	Optional : 200L/min	set	1
Carpet	Optional: Customize /Middle fleece/Light color	set	1



DR.HUGO

Your Oxygen Therapy Expert

POWER BY:



LANNX Bio&Med Co.,Ltd

Website: www.lannx.net/hbot/

E-Mail: sale@lannx.net

-Search our channel in YouTube to view operation video:

<https://www.youtube.com/@oxygensupplycenter>